

## Knee Suspension Sleeve

The **ISO Preferred Knee Suspension Sleeve** reduces the time it takes to heal discomfort from injuries or post-surgery. Please follow your orthopedic professional's instructions explicitly and adjust or discontinue use only as directed to achieve maximum benefit for your concerns.



ISO Preferred Knee Suspension Sleeve (KN208) is intended for use with hinged knee braces to provide additional support and comfort and to help prevent brace migration or skin irritation. The brace is effective in providing soft tissue compression to help aid the healing process.

1. The ISO Preferred Knee Sleeve may be applied in various ways to provide support, relief or comfort:

- a) Apply **UNDER** another brace. (1a). This application is intended to prevent slippage of the outer knee brace or to provide additional comfort or padding. Sleeve may be applied above the knee, at the knee or under the knee as needed or prescribed by your physician.

**While in a seated position**, apply sleeve directly to thigh, or directly over the knee or to the upper calf as prescribed by your physician. Pull the handle firmly for a snug fit and secure in place. **DO NOT OVERTIGHTEN.**

- b) Apply **OVER** another brace (1b). The sleeve may be applied over another knee brace to help secure or stabilize the brace.

**While in a seated position**, apply the sleeve by wrapping around the brace on upper, middle or lower section of the inner brace and securing for a snug fit. **DO NOT OVERTIGHTEN.**

2. The handle may be detached from the strap for cleaning. When reattaching, make sure the holes extend beyond the end of the sleeve and the labels are on the opposite side.


1a



1a



1b

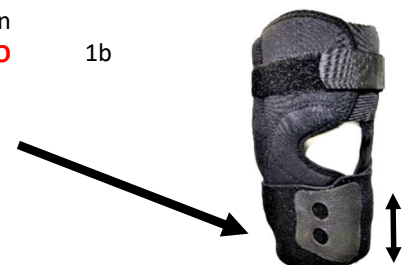

2





### INTRODUCTION

Your **ISO Preferred Knee Suspension Sleeve (KN208)** is designed to provide your body with the support and stability that it needs to successfully manage or eliminate your mobility concerns. It is critical that you follow your orthopedic professional's directions closely to ensure that the apparatus is worn properly and in the manner that will be most beneficial in addressing your issues. For correct application, please read the enclosed instructions and watch the on-line instructional video (if available) at [www.ISOpreferred.com](http://www.ISOpreferred.com) so that this device can provide you with superior long-term usability and care.

### USAGE RECOMMENDATIONS

It is vital that you read the included instructions and when applicable watch the instructional video for the **ISO Preferred Knee Suspension Sleeve (KN208)**. The brace is intended for application by qualified professionals only. It is not recommended for patients to fit or apply or adjust their new brace without an orthopedic professional's expertise, unless otherwise instructed. Further questions regarding adjustments or changes, frequency of use, and/or duration of use should be presented to your orthopedic professional, as they are familiar with the specific needs and conditions, for which this device was prescribed.

### CARE AND CLEANING INSTRUCTIONS

Please do not place your brace in the washing machine or dryer at any time. When preparing to clean your **ISO Preferred Knee Suspension Sleeve (KN208)** you will want to separate the parts where applicable and care for them individually.

- Do not dry clean, iron, bleach, or use an abrasive cleaning instrument on any part of the brace. Failure to follow these instructions may damage your brace and void the manufacturer's warranty and it may no longer work properly.
- Hand wash each section in cold water with a mild detergent. Lay flat to air dry.
- Allow all parts to air dry prior to reassembling the device.

### CAUTION

The **ISO Preferred Knee Suspension Sleeve (KN208)** is for single patient use only. While every effort has been made to obtain maximum strength, function and durability, there is no guarantee that injury will be prevented through the use of this device, especially if the instructions are not followed correctly. **DO NOT OVERTIGHTEN.**