

## Size Information

### Wrist & Forearm Splint



Wrist Sizing Guide-Measure above wrist bone.

#### The Maximus Wrist Brace

Item#	Size	Inches	CM
WB0104-148-01	Small	6 1/2 - 7 1/2	16-19
WB0104-148-02	Medium	7 1/2 - 8 1/2	19-21
WB0104-148-03	Large	8 1/2 - 9 1/2	21-24
WB0104-148-04	X-Large	9 1/2 - 10 1/2	24-27

## Washing Instructions

Hand wash in cold water using mild soap.  
Line dry only. Do not bleach and avoid direct  
heat or sun.



## Application Instructions

The Maximus Wrist Brace

### Caution


If you experience any pain, swelling, sensation changes, or any unusual reaction while using the product, please consult your physician immediately.




DR Medical  
203 Rio Circle Ste.B  
Decatur, GA 30030  
[www.drmedsupplies.com](http://www.drmedsupplies.com)

## INSTRUCTIONS FOR USE


### WB0104-148 The Maximus Wrist Brace

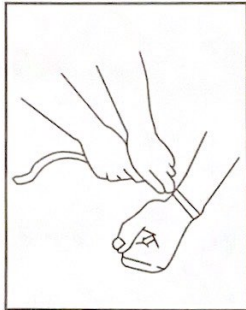


1. Use a soft / flexible tape measure.  
Note\* Measurements will be more accurate if you have someone measure for you. Place the tape measure around the wrist and measure the total circumference around your wrist. Use size chart to determine appropriate size brace based on wearers wrist measurements.

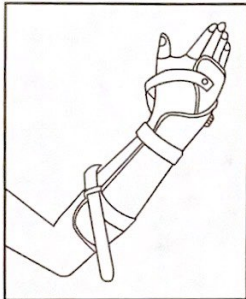

2. Loosen all straps of the brace for easier wrist and forearm insertion.

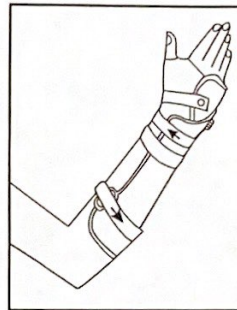

3. Pass thumb through thumb tunnel and refasten all straps to desired comfort and support level.

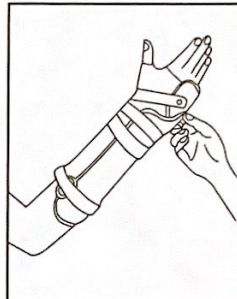

4. Adjust wrist hinge to desired wrist deviation. Once desired deviation is obtained, tighten screw to secure.